

# MURI2 Task 2 Transfer Protocol

Work is split into 10-day units. This is so we can guarantee that we have enough supplies for 10-day transfers.

**Goal: To ensure transparency in Task2**

## Shaker layout

There are two Innova 44 shakers in the Lennon Lab. Each has a maximum capacity of 96 50 mL Erlenmeyer flasks.

The top shaker lines that are transferred at a volume of 1 mL.

The bottom shaker contains lines of a "backup" species that are transferred at a volume of 1 mL.

Within the top shaker the lines are grouped into five blocks.

To briefly go over the layout, we have one transfer sizes, three transfer times, 7 strains (6 + 1 backup species), and replicates of each treatment. In the top shaker there are 90 lines, all of those lines are transferred at either the high or low dilution are contained within a shaker. The 90 lines are divided into 5 blocks (the number of replicates) and within each block the position of transfer time and strain permutations are randomized using a custom Python script that is on GitHub.

In the bottom shaker there are 25 lines, all of which are of the backup species.

In the backup shaker there are two replicates of each 100-day line of the 6 species, so 60 lines.

The layout of the top shaker is as follows:

|   | 1              | 2              | 3              | 4              | 5              | 6              | 7              | 8              | 9              | 10             | 11             | 12             |
|---|----------------|----------------|----------------|----------------|----------------|----------------|----------------|----------------|----------------|----------------|----------------|----------------|
| A | 0-T-A01-D1-J-1 | 2-T-A02-D1-C-1 | 2-T-A03-D1-A-1 | 0-T-A04-D1-B-1 | 1-T-A05-D1-A-1 | 0-T-A06-D1-C-1 | 2-T-A07-D1-A-2 | 0-T-A08-D1-D-2 | 1-T-A09-D1-D-2 | 1-T-A10-D1-C-2 | 1-T-A11-D1-A-2 | 1-T-A12-D1-J-2 |
| B | 1-T-B01-D1-C-1 | 1-T-B02-D1-F-1 | 0-T-B03-D1-F-1 | 1-T-B04-D1-D-1 | 0-T-B05-D1-A-1 | 0-T-B06-D1-D-1 | 2-T-B07-D1-J-2 | 1-T-B08-D1-B-2 | 0-T-B09-D1-F-2 | 1-T-B10-D1-F-2 | 0-T-B11-D1-B-2 | 0-T-B12-D1-J-2 |
| C | 2-T-C01-D1-J-1 | 2-T-C02-D1-D-1 | 1-T-C03-D1-J-1 | 1-T-C04-D1-B-1 | 2-T-C05-D1-B-1 | 2-T-C06-D1-F-1 | 2-T-C07-D1-F-2 | 0-T-C08-D1-C-2 | 2-T-C09-D1-D-2 | 2-T-C10-D1-C-2 | 2-T-C11-D1-B-2 | 0-T-C12-D1-A-2 |
| D | B1T            | B3T            | 1-T-D03-D1-A-5 | 2-T-D04-D1-J-5 | 1-T-D05-D1-B-5 | 0-T-D06-D1-D-5 | 0-T-D07-D1-F-5 | 2-T-D08-D1-F-5 | 1-T-D09-D1-F-5 | 1-T-D10-D1-C-5 | 0-T-D11-D1-C-5 | B5T            |
| E | B2T            | B4T            | 0-T-E03-D1-A-5 | 0-T-E04-D1-B-5 | 2-T-E05-D1-D-5 | 1-T-E06-D1-D-5 | 2-T-E07-D1-A-5 | 2-T-E04-D1-B-5 | 1-T-E09-D1-J-5 | 2-T-E10-D1-C-5 | 0-T-E11-D1-J-5 | B6T            |
| F | 0-T-F01-D1-J-3 | 1-T-F02-D1-C-3 | 1-T-F03-D1-F-3 | 1-T-F04-D1-D-3 | 0-T-F05-D1-B-3 | 0-T-F06-D1-A-3 | 2-T-F07-D1-F-4 | 1-T-F08-D1-C-4 | 0-T-F09-D1-F-4 | 1-T-F10-D1-B-4 | 2-T-F11-D1-A-4 | 0-T-F12-D1-B-4 |
| G | 2-T-G01-D1-B-3 | 1-T-G02-D1-A-3 | 1-T-G03-D1-J-3 | 2-T-G04-D1-A-3 | 0-T-G05-D1-C-3 | 2-T-G06-D1-F-3 | 2-T-G07-D1-J-4 | 2-T-G08-D1-C-4 | 0-T-G09-D1-C-4 | 2-T-G10-D1-B-4 | 1-T-G11-D1-A-4 | 0-T-G12-D1-A-4 |
| H | 2-T-H01-D1-J-3 | 0-T-H02-D1-D-3 | 1-T-H03-D1-B-3 | 0-T-H04-D1-F-3 | 2-T-H05-D1-C-3 | 2-T-H06-D1-D-3 | 1-T-H07-D1-D-4 | 0-T-H08-D1-D-4 | 1-T-H09-D1-F-4 | 1-T-H10-D1-J-4 | 0-T-H11-D1-J-4 | 2-T-H12-D1-D-4 |

The layout of the bottom shaker is as follows:

|   | 1              | 2              | 3              | 4 | 5 | 6 | 7 | 8              | 9              | 10             | 11             | 12             |
|---|----------------|----------------|----------------|---|---|---|---|----------------|----------------|----------------|----------------|----------------|
| A | 1-B-A01-D1-P-1 | 2-B-A02-D1-P-1 | 0-B-A03-D1-P-1 |   |   |   |   | 2-B-D02-D1-P-4 | 2-B-A02-D1-P-1 | 2-B-E02-D1-P-5 | 2-B-B03-D1-P-2 | 2-B-C02-D1-P-3 |
| B | 1-B-B01-D1-P-2 | 0-B-B02-D1-P-2 | 2-B-B03-D1-P-2 |   |   |   |   | 2-B-C02-D1-P-3 | 2-B-D02-D1-P-4 | 2-B-E02-D1-P-5 | 2-B-A02-D1-P-1 | 2-B-B03-D1-P-2 |
| C | 0-B-C01-D1-P-3 | 2-B-C02-D1-P-3 | 1-B-C03-D1-P-3 |   |   |   |   |                |                |                |                |                |
| D | 0-B-D01-D1-P-4 | 2-B-D02-D1-P-4 | 1-B-D03-D1-P-4 |   |   |   |   |                |                |                |                |                |
| E | 1-B-E01-D1-P-5 | 2-B-E02-D1-P-5 | 0-B-E03-D1-P-5 |   |   |   |   |                |                |                |                |                |
| F |                |                |                |   |   |   |   |                |                |                |                |                |
| G |                |                |                |   |   |   |   |                |                |                |                |                |
| H |                |                |                |   |   |   |   |                |                |                |                |                |

The layout of the backup shaker is as follows:

|   | 1              | 2              | 3              | 4              | 5              | 6              | 7              | 8              |
|---|----------------|----------------|----------------|----------------|----------------|----------------|----------------|----------------|
| A | 2-T-H01-D1-J-3 | 2-T-D04-D1-J-5 | 2-T-C07-D1-F-2 | 2-T-A07-D1-A-2 | 2-T-A03-D1-A-1 | 2-T-G07-D1-J-4 | 2-T-H06-D1-D-3 | 2-T-F07-D1-F-4 |
| B | 2-T-C10-D1-C-2 | 2-T-C11-D1-B-2 | 2-T-E07-D1-A-5 | 2-T-C01-D1-J-1 | 2-T-G04-D1-A-3 | 2-T-C05-D1-B-1 | 2-T-G06-D1-F-3 | 2-T-E05-D1-D-5 |
| C | 2-T-D08-D1-F-5 | 2-T-C02-D1-D-1 | 2-T-F11-D1-A-4 | 2-T-B07-D1-J-2 | 2-T-C09-D1-D-2 | 2-T-H12-D1-D-4 | 2-T-H05-D1-C-3 | 2-T-E10-D1-C-5 |
| D | 2-T-C06-D1-F-1 | 2-T-G10-D1-B-4 | 2-T-G01-D1-B-3 | 2-T-E08-D1-B-5 | 2-T-G08-D1-C-4 | 2-T-A02-D1-C-1 |                |                |
| E | 2-T-H05-D1-C-3 | 2-T-G06-D1-F-3 | 2-T-C10-D1-C-2 | 2-T-G01-D1-B-3 | 2-T-B07-D1-J-2 | 2-T-E05-D1-D-5 |                |                |
| F | 2-T-G04-D1-A-3 | 2-T-C07-D1-F-2 | 2-T-E10-D1-C-5 | 2-T-C01-D1-J-1 | 2-T-C09-D1-D-2 | 2-T-H01-D1-J-3 | 2-T-C11-D1-B-2 | 2-T-G08-D1-C-4 |
| G | 2-T-E07-D1-A-5 | 2-T-D08-D1-F-5 | 2-T-A02-D1-C-1 | 2-T-C06-D1-F-1 | 2-T-E08-D1-B-5 | 2-T-H12-D1-D-4 | 2-T-A07-D1-A-2 | 2-T-F11-D1-A-4 |
| H | 2-T-G07-D1-J-4 | 2-T-C02-D1-D-1 | 2-T-G10-D1-B-4 | 2-T-F07-D1-F-4 | 2-T-D04-D1-J-5 | 2-T-A03-D1-A-1 | 2-T-C05-D1-B-1 | 2-T-F11-D1-A-4 |

## Labelling system

Our labelling system differs from that used in Task 1, primarily because of the different experimental design. The primary goal in creating this label system was to place all relevant metadata in a label that is the same length across all variations of the label.

The label consists of 6 elements, 3 of which being numerical and 3 being alpha-numerical.

A table describing the order and element contents is below:

| Position | Describes         | Elements               |
|----------|-------------------|------------------------|
| 1        | Transfer time     | 0, 1, or 2             |
| 2        | Shaker            | T (top) or B(bottom)   |
| 3        | 96-well position  | A01 - H12              |
| 4        | Dilution          | D1                     |
| 5        | Strain            | A, B, C, D, F, J, or P |
| 6        | Replicate (block) | 1-5                    |

Below are tables describing what each element means in the labelling system that is not intuitively clear.

**Element 1: Transfer time.**

| Transfer time | Description       |
|---------------|-------------------|
| 0             | 1-Day Transfers   |
| 1             | 10-Day Transfers  |
| 2             | 100-Day Transfers |

**Element 4: Dilution.**

| Abbreviation | Description   |
|--------------|---------------|
| D1           | 1:10 dilution |

**Element 5: Strain**

| Abbreviation | Description                            |
|--------------|----------------------------------------|
| A            | <i>Arthrobacter</i> sp. KBS0703        |
| B            | <i>Bacillus subtilis</i> 168           |
| C            | <i>Caulobacter crescentus</i> NA1000   |
| D            | <i>Deinococcus radiodurans</i> BAA-816 |
| F            | <i>Pedobacter</i> sp. KBS0701          |
| J            | <i>Janthinobacterium</i> sp. KBS0711   |
| P            | <i>Psuedomonas</i> sp. KBS0710         |


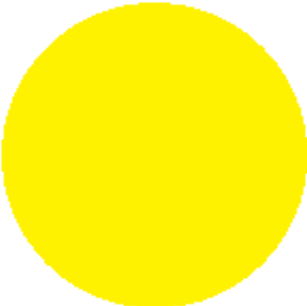
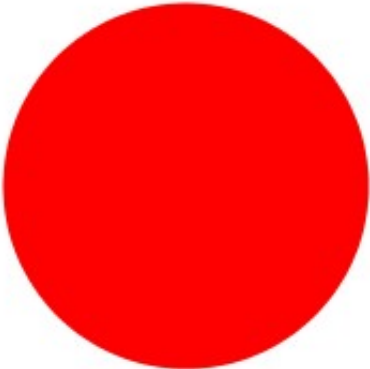
**Element 6: Replicates**

---

| Abbreviation | Description            |
|--------------|------------------------|
| 1            | Block 1 (top-left)     |
| 2            | Block 2 (top-right)    |
| 3            | Block 3 (bottom-left)  |
| 4            | Block 4 (bottom-right) |
| 5            | Block 5 (center)       |

**Transfer time**

The color of the sticker which the label is printed on indicates the transfer time of that line.

| 1-Day transfers                                                                    | 10-Day transfers                                                                   | 100-Day transfers                                                                   |
|------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|
|  |  |  |

**Example**

Say we have a line of Janthinobacterium that is being transferred every day. It is the fifth replicate and it has been randomly assigned to position G10 on the 96-well format shaker tray.

Our label would be **0-T-G10-D1-J-5**.

**Prep**

Every 10 days we will need a minimum of 360 50 mL Erlenmeyer flasks filled with 10mL of media. During 10-day increments where there is a 100-day transfer we need 420 flasks. It is good practice to make sure that there are at least ~ 20% more flasks with media at hand.

## Media

The recipe for the complex media created for Task 2 is PYE with 0.2% glucose and 0.1% casamino acids. The recipe is as follows:

- 2 g bactopectone
- 1 g yeast extract
- 0.3 g  $\text{MgSO}_4 \times 7 \text{ H}_2\text{O}$
- 2 g dextrose
- 1 g casamino acids

100  $\mu\text{L}$  of a 50  $\mu\text{g}/\text{mL}$  solution of cycloheximide is added to each flask to prevent fungal contamination

Bring up to 1 L final volume with DI  $\text{H}_2\text{O}$  Add 15 g Agar/ L if you want plates Sterilize by autoclaving

## Flasks

Labelled flasks for 10 days of transfers (5 2-day transfers and 1 10-day transfer) will be prepped every 10 days.

For prepping flasks, the protocol is as follows. Using a calibrated liquid dispenser, aliquot 10 mL of the MURI media into clean 50 mL Erlenmeyer flasks. Place a 30 mL plastic beaker over each flask. Then, sterilize a tray of those flasks marked with autoclave tape with  $\sim 3/4$  in. of water for 30 minutes in the autoclave. Once you've removed the flasks and they are cool, place them onto a dry tray and dump the water.

## Tubes

The same outline for flasks applies to tubes. Fill with 10 mL of media, autoclave with metal caps.

## Tips

We are using USA Scientific filtered pipet tips for our transfers. The Lennon Lab has a set of hood-dedicated pipettes, so those will always be available for transfers. Will will be ordering and keeping track of pipet tips over the course of the experiment.

## Printing labels

The IP address of the printer in the Lennon Lab is on the printer. To print the stickers, do the following:

- Choose sticker color by what transfer times you're printing (see above image). These stickers will be above the printer in the Lennon Lab.
- Open up the sticker format in the Box folder  $\sim/\text{MURI2}/\text{Aim2}/\text{Transfers}/\text{Labels}$ . To get all 60 labels for a

particular transfer time, you'll need to print three sheets of stickers.

- Open the slot on the printer marked with orange tape and insert your sticker sheet sticker side up.
- Print the stickers using the word document.

For a normal one day transfer you will need 35 stickers.

## **Setup**

Ten days of flasks will be labelled and organized before the start of a ten-day cycle. In each ten day cycle there are five one-day transfers and one ten-day transfer, so there a total of 6 trays that need to be prepped in a normal week (not including weeks where there are 100 day transfers). Each tray is labelled with a sticker at the position where the flask should be placed.

## **Washing**

Dump the contents of each used 50 mL flask into the large flask next to the sink in the Lennon Lab.

Fill each flask half-way or tube with water and autoclave to help remove any residual biofilms, then proceed with washing.

For washing the 360 plus 50 mL flasks that will be dirtied each week, we've purchased a 81 IXC 14 dishwasher rack from Lancer. This rack is compatible with the fourth floor dishwasher and is capable

The protocol for the dishwasher in on the fourth floor.

## **Transfers**

I have created a Google calendar detailing when transfers are occurring and it is available [here](#)

## **Transfer procedure**

When someone walks in to do a normal one-day transfer there should be a tray already prepared. Each tray is divided into two sections, each for one day. Within each section there are five rows of six. Each row corresponds to a block on the 96-well formatted tray and tells you what you are transferring that day. Double check to make sure that the flask order is correct.

Place your pipet tips and the vortexer in the hood and UV for ~10-15 minutes. Press the bottom on the hood that powers the outlet inside the hood.

Starting with the top shaker, remove the tray and then turn on the UV light. Place the lines you're transferring

that day on a new tray that has been wiped down with ethanol. Bring your tubes to the hood and ethanol your gloves.

Working across each block, transfer 1 mL of culture from the old flasks into the new ones.

Carry the tray of newly transferred lines over to the 96-well formatted tray and place each line into its 96-well position, indicated by the third object on the label on the sticker. Turn off the UV light and carefully place the tray back into the shaker.

Do the same thing for the bottom shaker.

Take all the old flasks and put them in a tray labelled with today's date. Take the flasks from the previous transfer and label them as waste. They'll be cleaned by (insert undergrad's name here).

## **Human error**

A goal with the protocol is to minimize the probability of mistakes occurring. However, if at any point you doubt the sterility of some supply during a transfer (ex. pipet tip brushes against glove, 10 mL dilution tube splashes out, etc), get a new one. There are always spares ready at hand.

If you notice that a flask is translucent (i.e. no growth), go back to the lines from either one or two days ago and transfer from that flask. Record the date, what line didn't grow, and how many steps back you had to transfer from (1 for one-days back, 2 for two-days back if the issue was with a one-day transfer).

## **Responsibilities**

Will Shoemaker

- Overseeing Task2
- Ordering
- Collecting samples
- generating and analyzing data

Evgeniya Polezhaeva

- Transfers
- Monitoring the lines for contamination
- Collecting samples

Bianca Blaettner

- Making media



- Prepping flasks
- Waste disposal
- Printing labels
- Washing glassware

## Emergency contacts

---

Will Shoemaker

- Phone: 703-554-9264
- E-mail: [wrshoema@uemail.iu.edu](mailto:wrshoema@uemail.iu.edu)

Dr. Jay T. Lennon

- [lennonj@indiana.edu](mailto:lennonj@indiana.edu)